



DS500 DIGITAL BATHROOM SCALE

# Instructions Manual Product Instructions

MODEL NUMBER: DS500  
QUICK INSTALL VIDEO



## Warning!

Do not overload the scale by exceeding the scale's capacity.

Do not use the scale with wet feet or in a wet environment.

Do not use it on any carpeted surface as that can affect the operation and the accuracy of the scale.

Only use the scale on a level, hard surface and ensure that the scale will not move when in use.

## Specification parameters

Weighing Capacity: 440x0.1 lb (200x0.1 kg)

Unit: lb, kg, st

Dimensions: 12x12x4

## Prompt

### Overload Condition




When the EEEE is displayed, the scale weight is past the maximum capacity of the scale. Only use the scale within the maximum range of the scale to prevent damage.



## Low Battery


The battery needs to be replaced.  
Replace the 2 AAA batteries.

# Operation and use

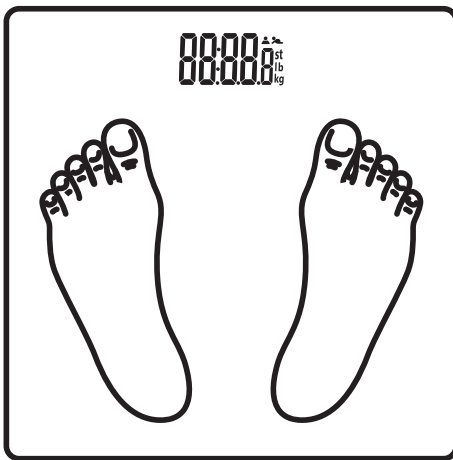


Step onto the scale. The scale will power on and provide the patient's weight. The scale will automatically turn off after 15 seconds while displaying the patient's weight.

# Unit conversion

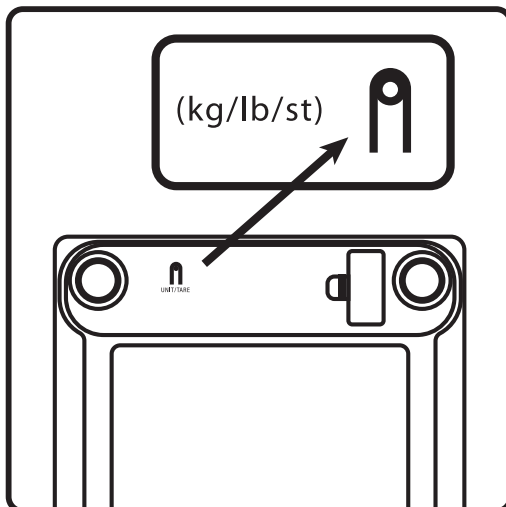


This product has three units available. Select the unit by turning over the scale and pressing the unit button to select the desired unit. Once selected this unit will be active until the selection is changed.

# Installing Batteries

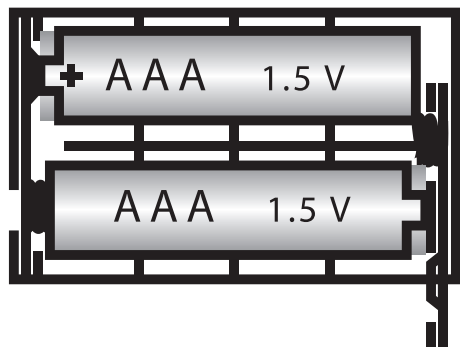

This product uses 2\*1.5V AAA batteries. These come installed, remove the shipping tag when installing the scale. When installing new batteries, be sure to install them with the correct polarity as indicated in the battery compartment.

## Mother Baby Mode

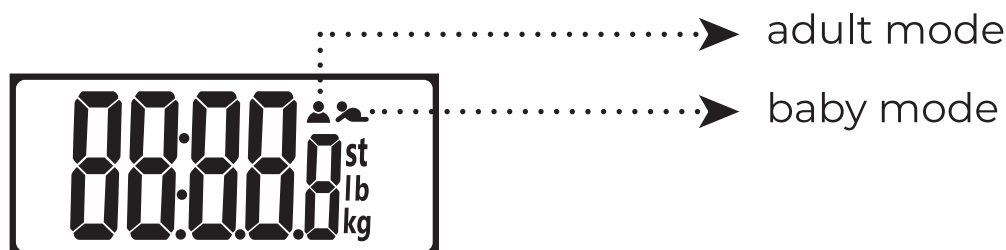

To place the scale in baby mode, turn on the scale by standing on the scale, while on press and hold the Unit selection button for three seconds. Once in baby mode, the scale will display the adult icon.

To exit the baby mode, turn on the scale by standing on the scale, while on press and hold the Unit Selection button for three seconds. The adult icon will no longer be displayed.